

S-CAPE® MEDICAL MULTICONSOLE D1

D1 · 24 | 27 | 32 | 42 | 47 | 55

other display sizes on request



Ready for the **PRIME365** integration solution.

Different surfaces and colors

The next generation MULTICONSOLE to manage your growing digital workflow needs.

The **MULTICONSOLE** is the centrepiece for your digital OR workflow.

Tailored to fit your needs with:

- optional display sizes and technologies
- a variety of input- and output add-ons
- a multitude of finishes and colors
- all contained in-wall or on-wall
- »zero footprint« deployment

Provides display of medical information and image data without compromising on the highest hygienic and ergonomic standards.

Further information about the **S-CAPE® MEDICAL MULTICONSOLE (MMC)** and the **PRIME365 software solution**:

S-CAPE® GmbH
Komturstraße 18a, 12099 Berlin
T: +49 30 7130297-0
E: info@s-cape.com
www.s-cape.com

technical data

console:

- + metal, powder coated RAL 9010 with non-reflecting safety glass for optimal protection of display surface
- + *optional*: RAL-colour of choice or stainless steel

displays:

- + full HD resolution (1920 x 1080 px), DICOM-preset from 32" to 55" display
- + *optional*: DICOM-preset for 24" and 27" display

powered by:

- + mainboard with INTEL® Core™ i5 processor, *optional*: INTEL® Core™ i7 processor
- + 8 GB RAM (expandable up to 32 GB)
- + LAN RJ 45 gigabit (10/100/1000mb/s),
- + *optional*: with galvanic isolation, glass fibre or WLAN
- + 250 GB SATA SSD *or other optional hard drives*

features:

- + solid swivel-mounted front folds upwards for service access
- + no visible cable routing
- + silicone short stroke-keyboard with touchpad, mouse tray and O.R. suitable silicone mouse
- + *optional*: completely disinfectable, antimicrobial coated, medical membrane keyboard
- + *optional*: with integrated CD/DVD drive*
- + USB 2.0 port front mounted (more USB ports optionally available)
- + ready for the **PRIME365** integration solution**

* loss of protection class IP65

** project-specific, additional housing for the technical equipment might be required



operated via touch and/or keyboard

+ Add-ons:


S-CAPE[®] MMC
EXTERNAL DISPLAY



CD/DVD-Drive*



Card-Reader*

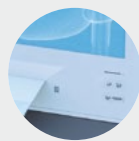


Multifunctional
Display



Loudspeaker



Barcode-/
QR-Code-Scanner

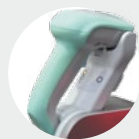

Interface
internal / external



The MULTICONSOLE is your entrance into the MMC family with PRIME365 software.

With the integration of additional modules you can enlarge your **MULTICONSOLE** to a central control unit in the O.R.**

Thanks to the new **PRIME365** software solution the surgical team always has the maximum amount of clinical information at hand. **PRIME365** arranges a more efficient flow of clinical data in and outside of the O.R. It is operated via touch and/or keyboard.

Your modular digital
O.R. integration solution.



- + Viewing of additional video sources
- + Video-routing
- + Image processing and archiving
- + Communication and/or integration with HIS-/PACS systems
- + Video-streaming and/or video-conferencing