


## S·CAPE® MEDICAL MULTICONSOLE D1-G

**D1-G** · 24 | 27 | 32 | 42 | 47

other display sizes on request

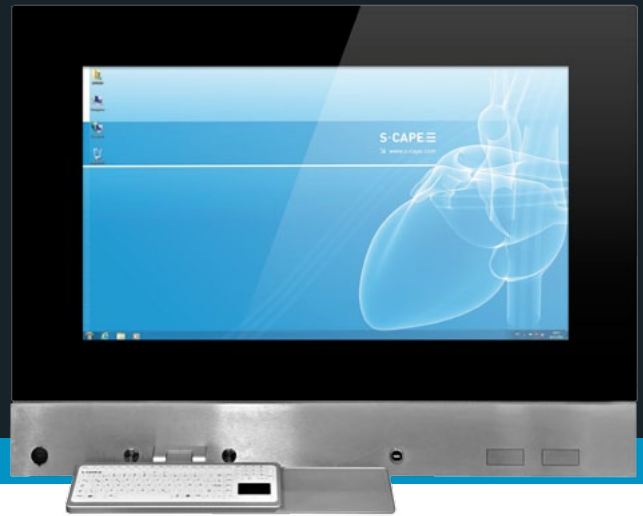

Ready for the **PRIME365** integration solution.

typical configuration

The new S·CAPE® glass front **MULTICONSOLE** features a one-piece glass front meeting the highest hygienic standards and optionally features multi-touch operation.

The glass front is resistant to cleaning agents and ageing. Contaminants are easy to remove from the smooth and hygienic front surface. The optional multi-touchscreen is useable with gloves, accepts 10-finger inputs simultaneously and 2-finger gestures.

### Tailored to fit your needs with:

- optional display sizes and technologies
- a variety of input and output add-ons
- flush installation
- self-contained in-wall enabling a »zero footprint« deployment

Provides display of medical information and image data.

**Further information about the S·CAPE® MEDICAL MULTICONSOLE (MMC) and the PRIME365 software solution:**

S·CAPE® GmbH  
Komturstraße 18a, 12099 Berlin  
T: +49 30 7130297-0  
E: info@s-cape.com  
www.s-cape.com

### technical data

#### console:

- + front surface: combination of non-reflective glass and stainless steel suitable for moist wipe disinfection
- + built-in base unit: metal, powder coated

#### displays:

- + full HD resolution (1920 x 1080 px), DICOM-preset from 32" to 47" display
- + *optional*: DICOM-preset for 24" and 27" display
- + *optional*: 24" or 42" multi-touch screen

#### powered by:

- + mainboard with INTEL® Core™ i5 processor, *optional*: INTEL® Core™ i7 processor
- + 8 GB RAM (expandable up to 32 GB)
- + LAN RJ 45 gigabit (10/100/1000mb/s), *optional*: with galvanic isolation, glass fibre or WLAN
- + 250 GB SATA SSD *or other optional hard drives*

#### features:

- + solid swivel-mounted front folds upwards for service access
- + no visible cable routing
- + silicone short stroke-keyboard with touchpad, mouse tray and O.R. suitable silicone mouse
- + *optional*: completely disinfectable, antimicrobial coated, medical membrane keyboard
- + *optional*: with integrated CD/DVD drive\*
- + USB 2.0 port front mounted (more USB ports optionally available)
- + ready for the **PRIME365** integration solution\*\*

\* loss of protection class IP65

\*\* project-specific, additional housing for the technical equipment might be required



operated via touch and/or keyboard

## + Add-ons:


S-CAPE<sup>®</sup> MMC  
EXTERNAL DISPLAY



CD/DVD-Drive\*



Card-Reader\*

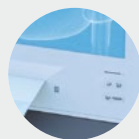

Multifunctional  
Display



Loudspeaker



Barcode-/  
QR-Code-Scanner



Interface  
internal / external



## The MULTICONSOLE is your entrance into the MMC family with PRIME365 software.

With the integration of additional modules you can enlarge your **MULTICONSOLE** to a central control unit in the O.R.\*\*

Thanks to the new **PRIME365** software solution the surgical team always has the maximum amount of clinical information at hand. **PRIME365** arranges a more efficient flow of clinical data in and outside of the O.R. It is operated via touch and/or keyboard.

Your modular digital  
O.R. integration solution.



- + Viewing of additional video sources
- + Video-routing
- + Image processing and archiving
- + Communication and/or integration with HIS-/PACS systems
- + Video-streaming and/or video-conferencing